



People.Health.Care.

Suprasorb[®] P sensitive Ag **NEW**

Gentle, versatile, protective
foam dressing with silver.



Suprasorb® P sensitive Ag **NEW**

Gentle, versatile, protective foam dressing with silver.

Suprasorb® P sensitive Ag is an absorbent foam dressing with antibacterial silver.

Key Benefits



Ionic Silver **extends wear time** by reducing bacterial (biofilm) colonization⁴



Prevents biofilm colonization for up to 7 days⁵



Dressing can remain on the wound site up to **7 days**⁵



Waterproof and **bacteria-proof** while allowing water vapor to escape⁵



Protective backing provides low coefficient of friction surface⁵



Can be worn **in the shower**⁵



Gentle yet secure silicone adhesive technology **minimizes patient pain** and trauma to the wound and skin during dressing changes^{1,6}

The new Suprasorb® P sensitive Ag foam dressing with antibacterial silver works with you and your patient to support your wound bed preparation, protecting against bacterial (biofilm) contamination, extending wear time and preventing medical adhesive related skin injury.

Indications

Suprasorb P sensitive Ag dressings are indicated for the management of moderate- to highly-exuding¹ acute or chronic wounds⁸, including:

- Venous leg ulcers
- Diabetic foot ulcers
- Pressure injuries
- Dehisced surgical wounds
- Graft donor and recipient sites
- Skin tears and abrasions

Suprasorb® P sensitive Ag **NEW**

Takes on the tough challenges of wound care.



Protection

- Protects from external moisture and microbial contamination
- Protects skin from medical adhesive-related skin injury
- Protects against biofilm buildup in the dressing to extend wear time




Efficient exudate management^{1,2}

- promotes quick absorption of exudate into the dressing
- maintains moist environment³





Softness & comfort³

- Gentle yet secure adhesion
 - Atraumatic — Pain-free removal
 - Reliable seal at the wound edge to protect skin
 - Extensible, flexible film and foam move with the patient for security and comfort.
- 

Suprasorb[®] P sensitive Ag **NEW**

The silver lining for sensitive wounds.

Suprasorb P sensitive Ag dressings with ionic silver provide you with the right combination of a moist wound environment¹ while managing moderate to highly exuding¹ wounds delivering effective antibacterial action to inhibit biofilm reconstruction.⁵

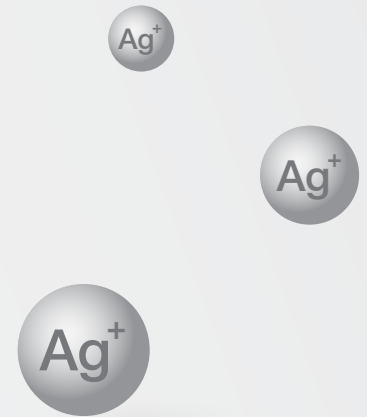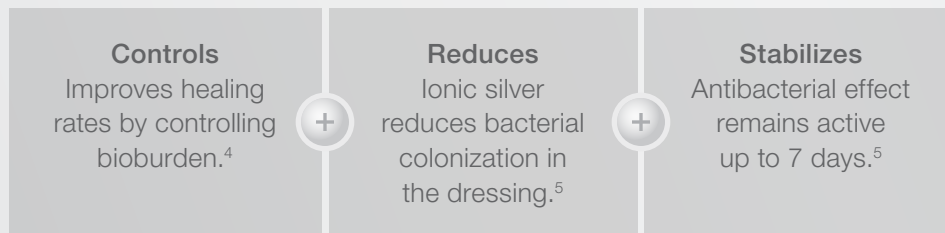

Suprasorb[®] P sensitive Ag

Five layers dedicated to wound healing.

Protective backing

- high moisture vapor transmission rate (MVTR)¹
- smooth and conformable³ with low coefficient of friction
- waterproof
- microbial barrier

Suprasorb P sensitive Ag border with super-absorbent core

- super-absorbent polymer fibers
- traps fluid under compression
- extends wear time and reduces risk of maceration²

Non-woven layer


- helps to distribute exudate

Absorbent polyurethane foam

- offers comfort and cushioning⁷
- efficient exudate management^{1,2}

OptiSil silicone wound contact layer

- perforations for swift exudate transfer
- soft and pliable for maximum patient comfort
- atraumatic release minimizes pain during dressing removal^{1,3}
- does not adhere to the wound³

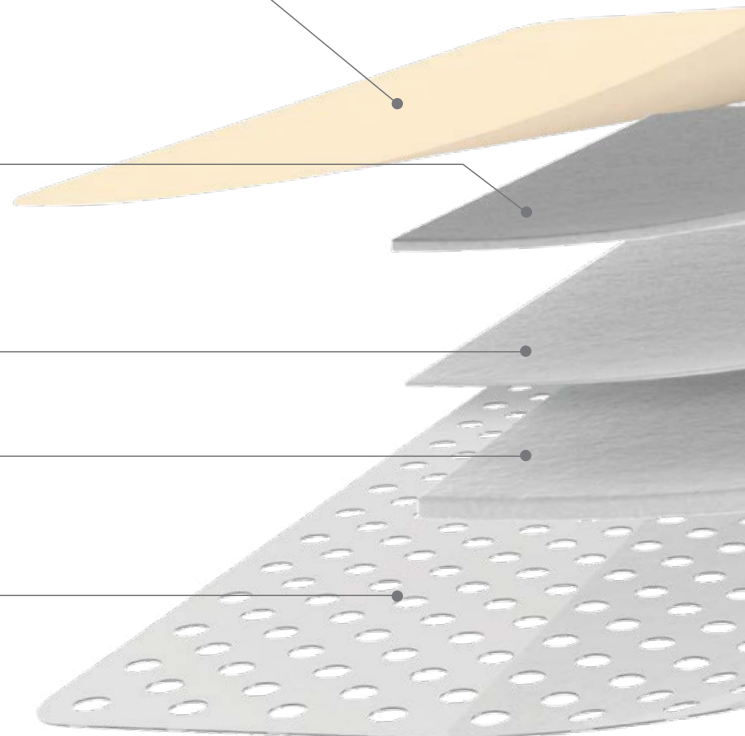

Suprasorb® P sensitive Ag NEW

Broad antibacterial action for all moderate to highly-exuding wounds.

border
with super-absorbent core



for moderately to heavily exuding wounds



non-border
without super-absorbent core



for moderately exuding wounds




Suprasorb® P sensitive Ag border

Indications


For the management of acute and chronic wounds such as⁸:

chronic wounds

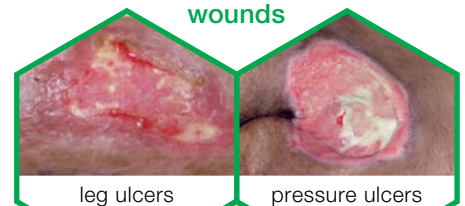

leg ulcers

pressure ulcers

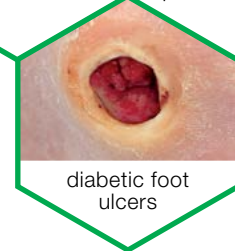

diabetic foot ulcers



skin abrasions



skin tears



wound dehiscence



skin graft donor & transplant sites



surgical staples / sutures



minimally invasive surgery

traumatic wounds

surgical wounds



People.Health.Care.



Suprasorb® P sensitive Ag border NEW
foam dressing with antibacterial silver and adhesive border
sterile, individually-sealed

Dressing Size	Pad Size	REF	HCPCS	Shipping Units (per box/case)
3 x 3 in. (7.5 x 7.5 cm)	1.8 x 1.8 in. (4.5 x 4.5 cm)	180010	A6212	10/200
4 x 4 in. (10 x 10 cm)	2.5 x 2.5 in. (6.2 x 6.2 cm)	180011	A6212	10/200
5 x 5 in. (12.5 x 12.5 cm)	3.4 x 3.4 in. (8.5 x 8.5 cm)	180012	A6212	10/200
8 x 8 in. (20 x 20 cm)	6 x 6 in. (15 x 15 cm)	180014	A6213	10/200

Suprasorb® P sensitive Ag non-border NEW
foam dressing with antibacterial silver
sterile, individually-sealed

Size	REF	HCPCS	Shipping Units (per box/case)
2 x 2 in. (5 x 5 cm)	180015	A6209	10/200
4 x 4 in. (10 x 10 cm)	180016	A6209	10/200
4 x 5 in. (10 x 12.5 cm)	180017	A6210	10/200
6 x 6 in. (15 x 15 cm)	180018	A6210	10/200
8 x 8 in. (20 x 20 cm)	180019	A6211	10/200

Sources

- Barrows, C. (2009). Enhancing Patient Outcomes- Reducing the Bottom Line: The Use of Antimicrobial Soft Silicone Foam Dressing in Home Health. Published in Home Healthcare Now, 27(5), 279-284.
- Pukki, T., Tikkanen, M., Halonen, S. (2010). Assessing Mepilex Boarder in post-operative wound care. Published in Wounds UK, 6(1), 30-40.
- Davis, P (2008). Evidence reviews: the clinical benefits of Safetac technology in wound care. Published in Journal of Wound Care.
- Barrett, S. (2009). Mepilex Ag: an antibacterial, absorbent foam dressing with Safetac technology. Published in November 2009 British Journal of Nursing (Mark Allen Publishing) 18(20):S28,S30-6
- Data on file.
- White, R. (2008). A multinational survey of the assessment of pain when removing dressings. Published in March 2008 Wounds UK 4(1)
- Stephannie K. Miller, Neal Sharma, Lauren C. Aberegg, Kimberly N. Blasiolo, Judith A. Fulton (2015). Analysis of the Pressure Distribution, Qualities of a Silicone Border Foam Dressing. J Wound Ostomy Continence Nurs. 2015;42(4):346-351. Published by Lippincott Williams & Wilkins
- Davies, P. McCarty, S. Hamberg, K. (2017), Silver-containing foam dressings with Safetac: a review of the scientific and clinical data, Published Journal of Wound, Vol. 26, No. Sup6a.

